

THE PLATFORM PARADOX: DEVICE DEVELOPMENT CHALLENGES WITHIN FAST EVOLVING DRUG DISCOVERY PIPELINES

Tim Quigg, Mariam Al-Amari and Joel Gresham of Crux Product Design discuss how to enable high-dose delivery for novel biologics, using *ex vivo* and *in silico* approaches to accelerate drug-device development and improve device performance, thus de-risking projects and saving costs.

There are a number of unmistakable factors that have propelled the success of platform devices. The ubiquitous push for at-home drug delivery has provided a clear and compelling business case for recalibrating drug discovery pipelines. Where biologic drugs were once almost exclusively delivered intravenously (IV), it is hard to ignore the increasing prevalence of subcutaneous (SC) assets across pharmaceutical pipelines. SC dosing regimens are now not only being developed in parallel to IV dosing, but increasingly being prioritised as the primary or standalone dosing option.

As many pharmaceutical companies normalise expansive pipelines with 100+ assets progressing through clinical milestones, the appeal of platform devices is clear; amortised device investments, with multiple assets benefitting from the same delivery device technology. While this “one-size-fits-all” elevator pitch is compelling, it may be in direct contrast to trends within large-dose biologics development.

Large-dose administration has seen the industry continue to challenge on-market

precedents, incrementally pushing the perceived limits of acceptable drug volumes, viscosities and delivery times. Furthermore, many large-dose assets exhibit highly complex rheological behaviours (Figure 1). At one extreme, ultra-high concentrations are breaching the 200 mg/mL precedents that have been established over the last decade, even reaching beyond 500 mg/mL. These trends have coincided with the development of non-aqueous suspensions, where a small modification in formulation chemistry can result in significant changes to syringeability, drug delivery and even pharmacokinetic (PK) performance.

Conversely, at the other edge of the design space, the industry is seeing comparably large doses being achieved with higher volumes (Figure 2). Whilst 2 mL autoinjector administrations have become the norm in recent years, pipelines are pushing far beyond this theoretical limit, with the advent of large-volume autoinjectors and on-body injectors proving to be a catalyst for this change. De-risking device development and

“DE-RISKING DEVICE DEVELOPMENT AND SELECTION WHILST WORKING WITHIN THIS EVER-CHANGING FORMULATION DESIGN SPACE HAS BECOME A KEY CHALLENGE FOR COMBINATION PRODUCT TEAMS.”

High Concentration	High Volume
Viscosity-reducing agents	High-volume (3 mL) autoinjectors
Ultra-concentrated suspensions	Hyaluronidase for >3 mL
Autoinjectors with powerful drive systems suitable for >20 cP viscosities	Injection pumps

Figure 1: Opportunities for high-dose delivery are commonly segmented into opposing formulation strategies – high concentration and high volume.¹

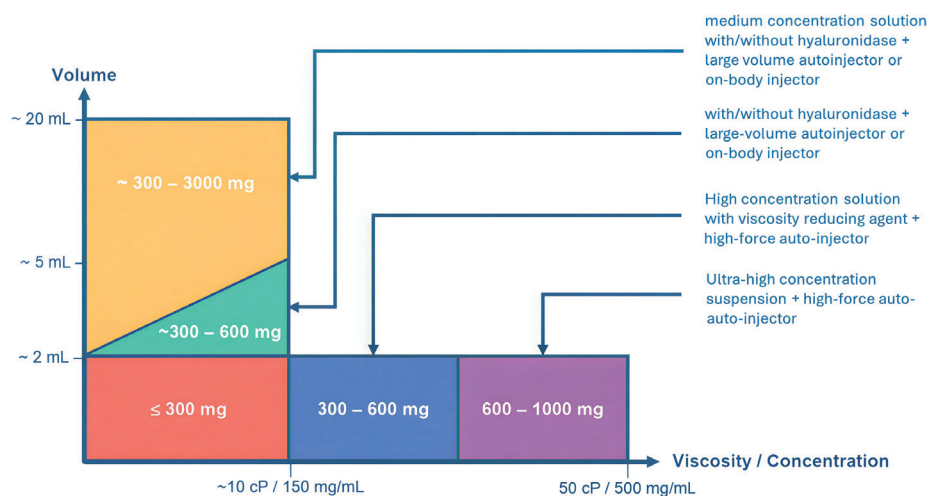


Figure 2: The high-dose design space is diverse, encompassing both higher volume and high concentration formulations, challenging a “one-device-fits-all” platform strategy.¹

selection whilst working within this ever-changing formulation design space has become a key challenge for combination product teams.

This can be described as “the platform paradox”. As off-the-shelf devices converge towards standardisation, next-generation biologics are pushing boundaries and allowing for increasingly unique, diverse design spaces. Thus, navigating this platform paradox requires novel research and development approaches, leveraging state-of-the-art analytical tools to make informed device decisions.

ENABLING HIGH-VOLUME DELIVERY

A variable that often remains unknown during device development is the behaviour of the injectate following injection. With the device landscape historically dominated by 1–2 mL or 1 cP prefilled syringes and spring-driven autoinjectors, it is no surprise that delivery profiles have often

been overlooked. A 1 mL injection yields negligible skin distention and therefore characterising bolus formation is rarely deemed essential for traditional SC delivery. However, when going beyond the traditional design space – expanding delivery volumes up to 20 mL and considering high-viscosity formulation technologies – the dynamic behaviours of SC bolus formation become significant.

Most biologics are absorbed into the vascular or lymphatic capillaries prior to ending up in systemic circulation. As there are multiple biological barriers to cross before reaching systemic circulation, and due to their size, some biologics show unfavourable bioavailability profiles after SC injection. For instance, some studies have shown that needle depth upon administration affects absolute bioavailability and PK parameters for biologics.² Therefore, this must be considered when developing injection devices for the SC route (Figure 3).

“CURRENT DEVICE PROGRAMMES CAN BENEFIT FROM SUPPLEMENTARY EX VIVO TISSUE TESTING, FRONT-LOADING THE CURRENT UNDERSTANDING OF DRUG-DEVICE-TISSUE INTERACTIONS PRIOR TO PRECLINICAL AND CLINICAL STUDIES.”

As bolus formation and PK performance become salient questions, traditional evaluation methods can prove both limited and costly. Preclinical and PK studies may be inevitable as a pipeline asset progresses towards market, but significant value can be derived through the adoption of *ex vivo* testing and advanced imaging. Where historic development programmes relied heavily on injections into air to demonstrate credible delivery times, current device programmes can benefit from supplementary *ex vivo* tissue testing, front-loading the current understanding of drug-device-tissue interactions prior to preclinical and clinical studies.

Despite traditional verification testing standards not requiring *ex vivo* evaluations, it is crucial to visualise depot formation and drug dispersion post injection. This, alongside measurement of tissue back-pressure, can help to guide device design (e.g. injection rate, needle selection, injection depth, push-on-skin pressure) and de-risk projects progressing towards the clinic.

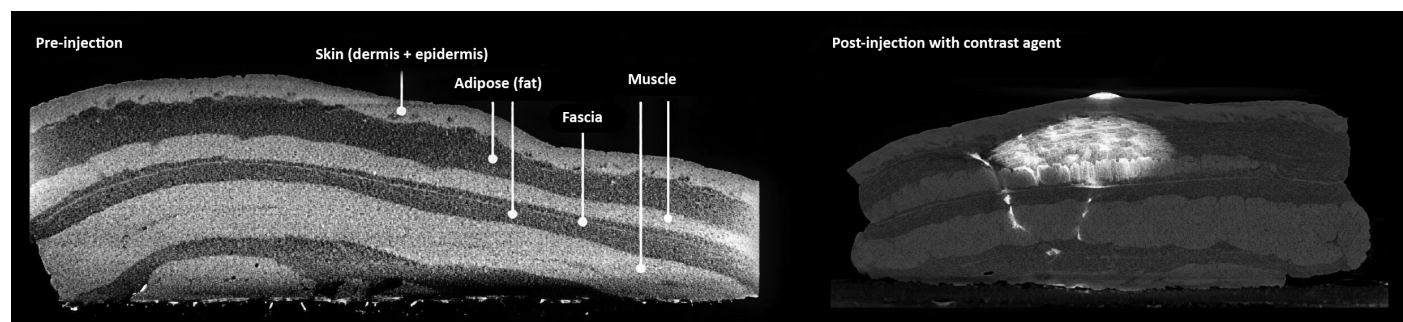


Figure 3: The SC region consists of adipose fat tissue, sometimes separated from muscle by the presence of fascia.

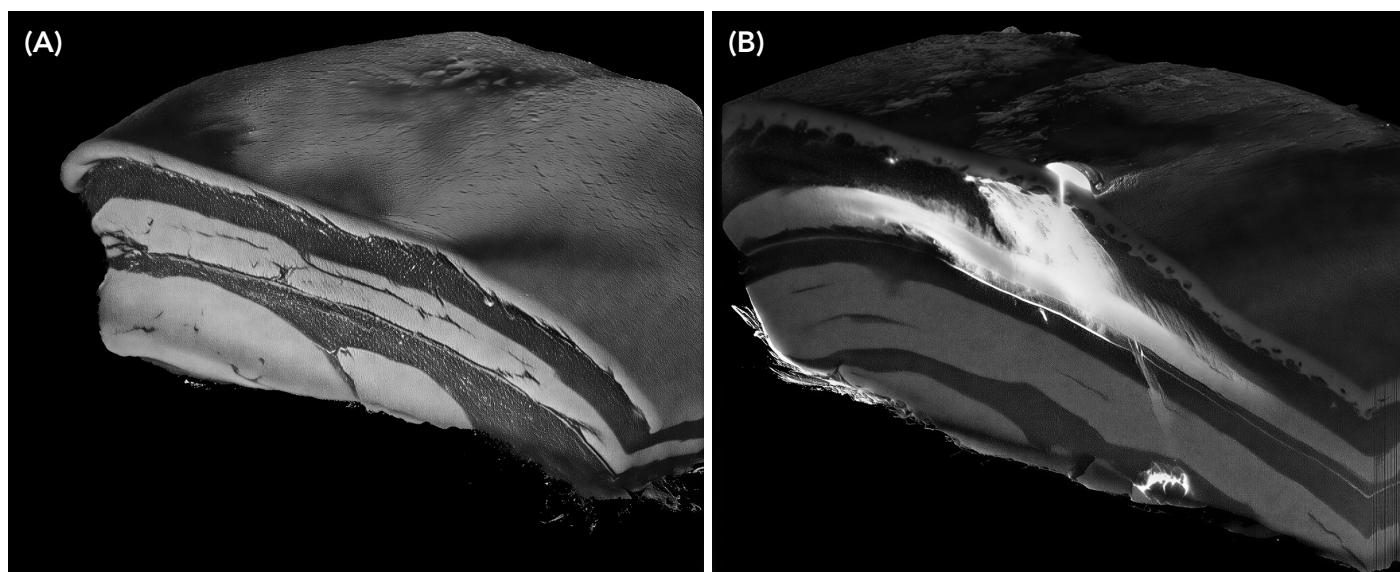


Figure 4: 5 mL bolus in *ex vivo* tissue, showing partial intramuscular penetration below the intended SC layer and significant leakage from the injection site. (A) Before injection. (B) After injection.

Public domain literature presents the results of *ex vivo* testing on porcine abdominal tissue, leveraging micro-CT scanning technology.³ The porcine model was selected for its translational relevance to human tissue, with resemblant skin thickness and skin tissue architecture. While *ex vivo* tissue does not account for full injectate absorption due to its lack of lymphatic drainage and wider physiological behaviours, it is recognised that minimal fluid removal occurs within the first 20 minutes post injection, enabling an accurate study of injectate dispersion in SC tissue immediately after injection.³

A 3D injected bolus volume can be derived following injection, providing datasets that can be interrogated to identify statistically significant variables (Figure 4). Data from these 3D images can be used to assess industry-wide areas of uncertainty, from predicting bioavailability changes as delivered volumes increase, to how bolus morphology may influence patient tolerability and pain.

Early-stage studies began as qualitative, comparative evaluations of in-development formulations. Today, this approach has matured, with refined test protocols, tissue handling frameworks, imaging procedures and data processing for investigating bolus spread, injection site reactions, leakage and skin distention. With advancing research, it has been found that this data can feed directly into device development programmes to:

- Inform optimal injection sites for different patient populations, which has an immediate effect on drug PK
- Study pain and tolerability as a function of bleb formation and size to design injection devices that may improve patient adherence
- Evaluate how injection rate and viscosity can influence leakage after needle removal, resulting in injectate loss, which influences how long the injection duration should be
- Assess device-skin interactions and how large-volume administration may affect tissue backpressure.

Combining the data derived from micro-CT imaging to visualise the injectate dispersion post injection has proven to be a valuable method to enhance due-diligence processes and inform the design space for in-development combination products.

MODELLING THE HIGH-CONCENTRATION DESIGN SPACE

In silico modelling and simulation technologies can be used to inform and establish the design space for high-dose SC biotherapeutics. The high-concentration domain remains poorly defined due to the novelty of high-dose biotherapeutic formulations and a historical focus on small molecules. These typically do not exceed 2 mL dose sizes and remain below a concentration of 250 mg/mL. As formulation technologies mature, so too must analytical techniques.

Advanced simulation approaches using finite element analysis and computational fluid dynamics (CFD) can be deployed to create a digital twin to simulate the complex mechanics, material behaviours and fluid interactions of high-concentration therapies. Digital twins are already acknowledged for their ability to predict complex device functionality, though their utility may be extended to include drug-device interactions. This enables device engineers to collaborate

"DATA FROM THESE 3D IMAGES CAN BE USED TO ASSESS INDUSTRY-WIDE AREAS OF UNCERTAINTY, FROM PREDICTING BIOAVAILABILITY CHANGES AS DELIVERED VOLUMES INCREASE, TO HOW BOLUS MORPHOLOGY MAY INFLUENCE PATIENT TOLERABILITY AND PAIN."

closely with formulation teams, evaluating how in-development drug assets may perform with available platform devices. Further challenging the platform paradox, these digital approaches can be leveraged to define clear device performance limits across the emerging formulation design space.

Amongst today's high-dose formulation technologies, suspensions remain prevalent. Whilst overall viscosities may be kept relatively low, suspension-based formulations are not without complexity when integrating them with platform devices. The most pertinent risk relates to poor syringeability (drawing into a syringe from a vial) and injectability (delivering out of a syringe into tissue) due to particle-induced needle occlusion. These issues stem largely from inter-particle interactions, including sedimentation and agglomeration. Whilst a fully mixed suspension avoids needle clogging, the risk is increased when a formulation is not fully reconstituted or resuspended before use (Figure 5).

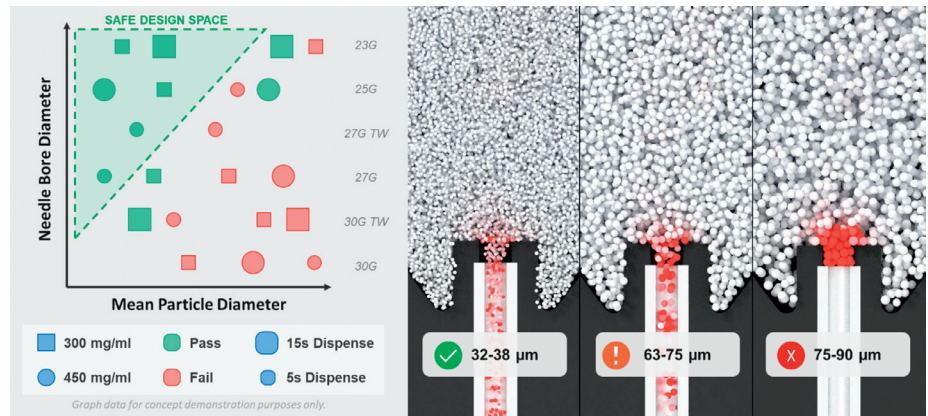


Figure 5: CFD and discrete element modelling to optimise drug particle and needle size specifications.

Addressing the significant complexity of biologics, with a diverse range of formulation variables, by conducting a comprehensive design-space investigation into needle occlusion would be costly and time-consuming. In contrast, a large virtual design-space study can be completed in a matter of weeks by leveraging computational models

and high-performance computing. Hundreds of simulations can be run in parallel, exploring thousands of device-formulation combinations to inform development strategies. Such models can evaluate how factors such as particle morphology, cohesivity and vehicle rheology can directly influence needle occlusion. Combining these studies with

GORE® IMPROJECT®

Silicone-Free, 0.5 mL Plungers for Bare (non-siliconized) Prefilled Syringes

Injection innovation begins with **BLUE!**



NEW PREFILLED SYRINGE PLUNGER SOLUTION FOR DELIVERY OF SILICONE-FREE* INTRAVITREAL INJECTIONS

Protect Complex, Sensitive Biologics and Gene Therapies in Prefilled Syringes

Gore focused on addressing complications with siliconized intravitreal prefilled syringe plungers:

- Requirements with significantly lower particle counts driven by <USP>789
- Silicone oil and floaters
- Silicone-induced inflammation
- Intraocular pressure (IOP) increases with repeated injections
- Interactions between silicone & polysorbates (fatty acids & particles)

* Silicone is not a raw material and is not added during the manufacture of the plungers

Together, improving life

pharmbio@wlgore.com

gore.com/pharmbio



the variability observed across the device landscape (e.g. needle size and inlet geometry) provides unparalleled insight into drug-device interactions when changing formulation properties.

ABOUT THE COMPANY

Crux Product Design is a technology-driven consultancy specialising in drug delivery and healthcare innovation. The company partners with leading pharmaceutical companies to develop inhalers, injection devices, wearables and novel delivery systems. Crux's multidisciplinary team combines engineering, sciences, human factors and design expertise to turn complex challenges into successful products. Certified to ISO 13485, ISO 9001, and ISO 14001, with an Ecovadis Gold rating, Crux is trusted by the world's top pharma companies. Based in Bristol, UK and Cambridge, MA, US, Crux delivers evidence-based innovation that drives the future of drug delivery.

REFERENCES

1. Bruin G et al, "An industry perspective on clinical development and regulatory strategies for subcutaneously administered high-dose biologics". *J Control Release*, 2025, Vol 386, art 114156.
2. Hu P et al, "Systematic review of device parameters and design of studies bridging biologic-device combination products using prefilled syringes and autoinjectors". *AAPS J*, 2020, Vol 22(2), art 52.
3. Gresham J et al, "Visualisation and quantification of subcutaneous injections of different volumes, viscosities and injection rates: An ex-vivo micro-CT study". *J Pharm Sci*, 2024, Vol 113(12), pp 3447–3456.

"HUNDREDS OF SIMULATIONS CAN BE RUN IN PARALLEL, EXPLORING THOUSANDS OF DEVICE-FORMULATION COMBINATIONS TO INFORM DEVELOPMENT STRATEGIES."



Tim Quigg

Tim Quigg, Development Director at Crux Product Design, supports development strategy and novel problem solving across drug delivery device development programmes. Crux collaborates with 18 of the top 20 pharmaceutical companies, providing Mr Quigg with foresight into cross-industry challenges. With a master's in Mechanical Engineering and a career spanning both combination product and medical device development, Mr Quigg champions front-loading of technical understanding, leveraging emerging tools to de-risk design, development and commercialisation.

T: +44 117 300 9788

E: tim.quigg@cruxproductdesign.com



Mariam Al-Amari

Mariam Al-Amari, Sales Executive at Crux Product Design, holds an MSc in Model-Based Drug Development and a BSc in Pharmacology from the University of Manchester (UK). At Crux, she applies her scientific expertise and industry insight to help pharmaceutical and medical device clients navigate innovation and development challenges, ensuring that solutions are informed, efficient and aligned with real-world needs.

T: +44 117 300 9788

E: mariam.al-amari@cruxproductdesign.com



Joel Gresham

Joel Gresham is Chief Scientist at Crux Product Design, leading a multidisciplinary team of technical consultants to deliver high-impact solutions to the world's largest pharmaceutical companies. Mr Gresham has led a range of medical technology and drug delivery projects spanning device and formulation design for companies including GSK, Novartis, Pfizer and Gilead.

T: +44 117 300 9788

E: joel.gresham@cruxproductdesign.com

Crux Product Design

Flatiron Building, 332-336 Paintworks, Bristol, BS4 3AR, United Kingdom
www.cruxproductdesign.com

BRAND NEW TITLE
LAUNCHING MARCH 2026

